


# EverPatch+

Permanent Scleral Reinforcement Matrix

2

Configurations  
in 1 pack

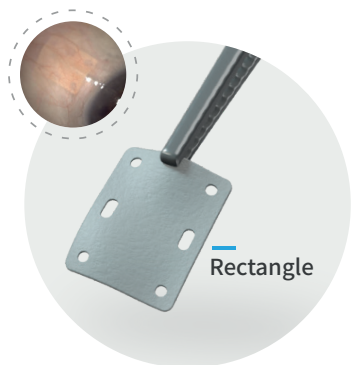

EverPatch+ is a Permanent Scleral Reinforcement Matrix - the first inert, synthetic, and non-degradable tissue substitute for ocular surface surgeries. It eliminates the need for donor tissues like cornea, sclera, and pericardium, which can degrade over time and carry disease risks.


Sterile and made from a non-woven polymer matrix, EverPatch+ integrates seamlessly with surrounding tissue by triggering cellular proliferation. It is indicated for scleral reinforcement and ocular surface reconstruction.

The EverPatch+ package includes two configurations:

**Rectangular shape:** For general scleral reinforcement and ocular surface reconstruction.

**Shield shape:** Optimized for reinforcing the ocular surface during glaucoma drainage device implantation, designed to reduce the patch's footprint at the limbus.

EverPatch+ meets the criteria of being permanent, sterile, immunologically inactive, cosmetically acceptable, and readily available.

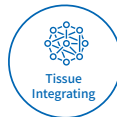

Order now and elevate your practice!

## Tapered Anterior Edge

Reduces the anterior footprint of the patch, maintaining distance from the edge of fornix-based flaps

## Anterior Stopper Band

Prevents anterior migration of the patch, maintaining distance from the edge of fornix-based flaps

## 100µm thickness

Conforms to the ocular surface anatomy thus reducing surface irregularities

## Bio-Stitching/Suturing Holes

Affix the device to the surgical site between sclera and conjunctiva, and enable accurate and simple suturing without the risk of cheesewiring

## Posterior Fixation Band

Stabilizes the patch on a glaucoma tube shunt – prevents lateral migration and reduces OR time and sutures required



For Inquiries and orders - [orders@corneat.com](mailto:orders@corneat.com)



“Thinner than a scleral patch, flexible, durable and optimized for manipulation & suturing. It's almost transparent under the conjunctiva with a more natural appearance” / **Dr. Ike Ahmed**





“Easy to implant and place sutures. Hydrated well” / **Dr. Jason Bacharach**

“Thinness an advantage” / **Dr. Glenn Graves**


# A NEW ERA OF MEDICAL DEVICES

## ABOUT US

CorNeat Vision is a leading and innovative ophthalmic medical device company, with a product portfolio addressing unmet needs in cornea, glaucoma, and ocular surface surgery. Our products, which include an artificial cornea, a synthetic scleral patch, and a revolutionary glaucoma drainage device, leverage our novel and proprietary permanent tissue-integrating material technology, the EverMatrix™, to seamlessly integrate with surrounding living tissue, offering enhanced, long-lasting solutions that reverse and prevent blindness.

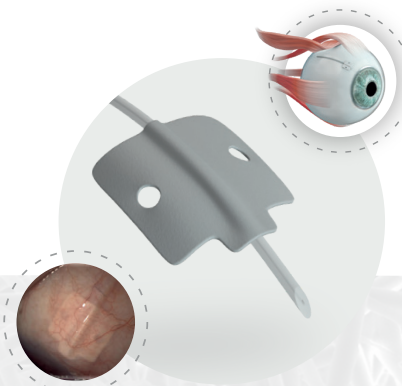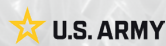

**eShunt**  
Drains Intraconally

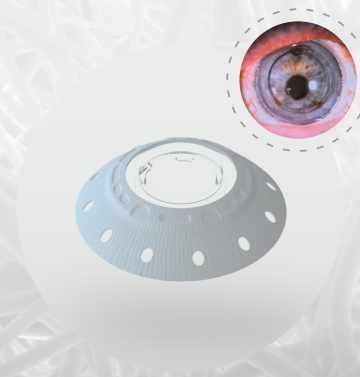

**CorNeat KPro**  
Revolutionary Artificial Cornea



Not available for sale